



# AD226 Diagnostic Audiometer

---

The portable true  
hybrid audiometer



# Stationary or portable – you decide!



---

**The AD226 is a powerful and fully portable hybrid audiometer.**

---

## **A portable true hybrid audiometer**

With the AD226 you can easily perform routine air and bone conduction threshold audiometry, but it also comes with multiple advanced features for other applications.

With a weight of only 1.3kg/2.9lbs. and compact footprint of only 30x23x9cm, 12x9x4 inches. it is ideal for portable application as well as small office places.

---

## **Familiarity + functional design**

The ergonomic design of the AD226 is exceptional with a high-resolution 4.3 inch color display that presents crystal clear images for easy to read audiograms. The familiar keypad design enables you to acclimate quickly without losing efficiency or speed.



---

## **A versatile True Hybrid audiometer**

The AD226 is a fully self-contained standalone audiometer that can run independently from a computer. It is also a fully PC-integrated audiometer (with the optional Diagnostic Suite software) offering all the benefits of a PC-based audiometer with full Noah, OtoAccess™ and EMR/HIS integration.

The AD226 is small for portability, yet presents a high professional image. The power of the internal processor provides PC-like functionality without the PC. Upload multiple patient data into the AD226 before taking it on the road, then store test results to each specific client and download the information to the database program of your choice (NOAH, OtoAccess™) when you get back to the office.



---


### Two versions: Basic and Extended with special test capabilities

The AD226 basic version (License B) includes air and bone conduction, masking, Stenger, Hughson Westlake (auto testing) and the ABLB test. In addition to this the AD226 extended version (License E) includes the Békésy, SISI, and Langenbeck (Tone in noise).

---

### A preview of the benefits

- Robust, fast and easy to use
- Built-in graphical display
- Storage for 500 patients
- Only 1.3kg / 2.9lbs.
- Direct print to PC-printers
- Possibility to run on batteries
- Optional PC integration (hybrid mode) with OtoAccess™, Noah storage and EMR interface capabilities

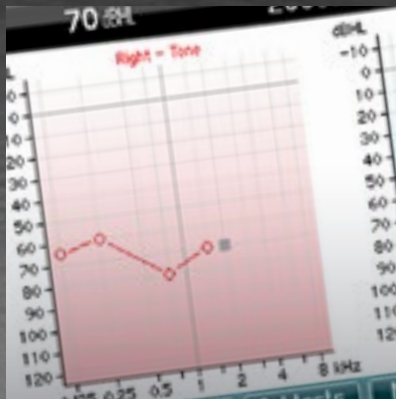

AD226 tone audiogram

- 
- Air and bone conduction audiometry with pure tone, warble and pulsed tones, WN and NB masking
  - Built-in talk forward microphone





---

## Interacoustics a/s

Interacoustics is a world leading diagnostic solutions provider in the field of hearing and balance assessment. We help the professional audiologic world reach new milestones through continuous developments and a constant focus on integration and direct customer value.

---

## With you at all times

Whenever you purchase audiological equipment with the Interacoustics brand on it, you are not just guaranteed a quality product but also highly trained support personnel with specialized knowledge. We operate in over 100 countries worldwide through a carefully selected network of distributors and service centres to ensure that you receive total support and service.

---

## Product specifications

All technical and hardware specifications concerning the AD226 can be downloaded from our website.

Read more at  
[interacoustics.com](http://interacoustics.com)

---

## Interacoustics a/s

Drejervaenget 8  
DK-5610 Assens

T +45 6371 3555  
F +45 6371 3522

info@interacoustics.com  
interacoustics.com

For sales and service in your area:



**Glenn Johnson**  
0414 452 007  
glenn@hearingcompany.com.au

Go online to  
explore our  
full product  
range

