

# Ezi Booth

Compact and portable

The Ezi Booth is our  
entry level, single person  
audiometry booth.


## Key Features

- Compact, single person booth
- Very good acoustic performance
- Portable with heavy duty castors
- Over 100 colour options
- Delivered fully assembled, Australia wide
- Designed and manufactured in Australia

**See next page for more details.**



# Compact and portable Ezi Booth

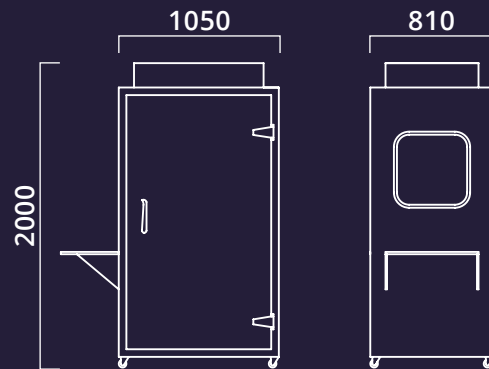

## High noise reduction

- Double-glazed window, high quality acoustic seals, double-wall steel construction
- Wall cavities packed with high density, sound-absorbent acoustic padding
- Designed to meet Australian Standard AS1269.4 specification
- Factory tested for internal noise levels prior to dispatch
- Perforated steel interior reduces reverberation (echo)

## Quality construction

- Solid construction with high quality Australian steel
- Durable paint, components and fittings
- Mounted on heavy duty castors
- Manufactured to suit a range of locations (mine sites, medical centres, research labs, etc.)
- Designed and built by Australian acoustic fabricators with over 30 years' experience

## Portable

- Arrives fully built on custom skid/palette, ready to unpack and move into your practice
- Easily move between rooms within a practice
- Two strong people can push over hard, flat surfaces
- Lift with a forklift/hoist onto a truck and relocate to another building



Cable connection panel



Internal equipment hooks

## Design features

- Standard 6.35mm sockets suit most audiometers
- External shelf (add or remove as required)
- Shelf slot for easy connection and less cable clutter
- Choose left or right door opening to suit room layout
- Wide selection of colours in hard-wearing finish
- Fits through standard doorway
- Wired to Australian Standard specification, fitted with circuit breaker and safety switch

## Which booth is right for you?

	Ezi Booth (most popular)	Flexi Booth	Maxi Booth
External dimensions	Width: 1050mm Depth: 810mm Height: 2000mm	Width: 1050mm Depth: 980mm Height: 2100mm	Width: 1100mm Depth: 1100mm Height: 2280mm
Internal dimensions	Width: 890mm Depth: 640mm Height: 1670mm	Width: 950mm Depth: 840mm Height: 1910mm	Width: 900mm Depth: 900mm Height: 1940mm
Wall thickness	50mm	50mm	100mm
Weight	350kg	250kg	500kg
Assembly	Delivered fully assembled	3 hours, 2-3 people	4 hours, 3-4 people
Portable	On castors - can be wheeled between rooms or buildings	Disassembly required for relocation	Disassembly required for relocation

Power	240 volt 10 amp (mains power), circuit breaker, safety switch
Door opening	Choose left or right opening
Inclusions	5m extension cable, external shelf, internal fan, low-heat LED light, carpet, equipment hooks
Options	Extended warranty, maintenance kit, USB socket, additional window in door

For sales and service in your area:

Glenn Johnson

☎ 1300 118 222  
✉ [contact@hearingcompany.com.au](mailto:contact@hearingcompany.com.au)