

Flexi Booth

Professional and spacious

The Flexi Booth is our
mid-range, single person
audiometry booth.


Key Features

- Larger internal space for comfort
- Medium weight hybrid panel construction
- Acoustic fabric interior
- Flat packed, ready for straightforward assembly
- Carry panels upstairs or around tight corners
- Australian designed and built

See next page for more details.



Professional and spacious Flexi Booth

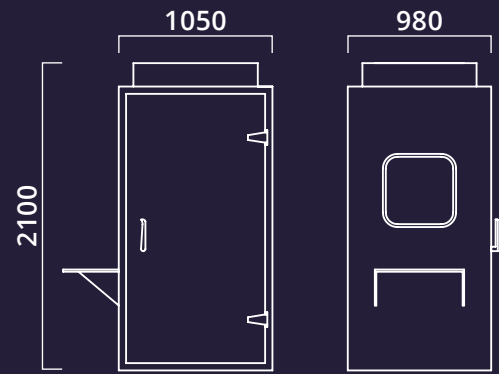

Excellent noise reduction

- Acoustic fabric interior for high noise absorption
- Double-glazed window, quality acoustic seals, double-wall steel construction
- High density, sound-absorbing acoustic panels
- Designed to meet Australian Standard AS1269.4 specification
- Factory tested for internal noise levels prior to dispatch

Quality construction

- Solid construction using quality Australian steel
- Commercial grade interior acoustic fabric
- High-quality components, fittings, seals and paint
- Manufactured to suit a range of locations (mine sites, medical centres, research labs, etc.)
- Designed and built by Australian acoustic fabricators with over 30 years' experience

Flexible and comfortable

- Medium weight hybrid steel/fabric panels can be carried up stairs or around tight corners
- Mid-level interior space offers more knee and leg room
- Fabric interior creates a softer, "corporate" look
- Easy entry/exit: low step, wide door opening, high ceiling
- Very stable, floor-mounted structure
- Optional door window for added light and visibility



Acoustic fabric interior



External equipment shelf

Design features

- Flat-packed tongue and groove design, for straightforward onsite assembly
- Can be disassembled, relocated and reassembled
- Standard 6.35mm sockets suit most audiometers
- External shelf included (add remove as required)
- Shelf slot for easy connection and less cable clutter
- Choose left or right door opening to suit room layout
- Wide selection of colours in hard-wearing finish
- Wired to Australian Standards specification with circuit breaker and safety switch

Which booth is right for you?

	Ezi Booth (most popular)	Flexi Booth	Maxi Booth
External dimensions	Width: 1050mm Depth: 810mm Height: 2000mm	Width: 1050mm Depth: 980mm Height: 2100mm	Width: 1100mm Depth: 1100mm Height: 2280mm
Internal dimensions	Width: 890mm Depth: 640mm Height: 1670mm	Width: 950mm Depth: 840mm Height: 1910mm	Width: 900mm Depth: 900mm Height: 1940mm
Wall thickness	50mm	50mm	100mm
Weight	350kg	250kg	500kg
Assembly	Delivered fully assembled	3 hours, 2-3 people	4 hours, 3-4 people
Portable	On castors - can be wheeled between rooms or buildings	Disassembly required for relocation	Disassembly required for relocation

Power	240 volt 10 amp (mains power), circuit breaker, safety switch
Door opening	Choose left or right opening
Inclusions	5m extension cable, external shelf, internal fan, low-heat LED light, carpet, equipment hooks
Options	Extended warranty, maintenance kit, USB socket, additional window in door

For sales and service in your area:

Glenn Johnson

☎ 1300 118 222

✉ contact@hearingcompany.com.au