

Maxi Booth

High performance

The Maxi Booth is our largest,
highest performing single
person audiometry booth.


Key Features

- Spacious interior
- Extra thick walls for high level noise reduction
- Sturdy floor mounted structure
- Flat packed, ready for straightforward assembly
- Over 100 colour options
- Delivered Australia wide

See next page for more details.



High performance Maxi Booth

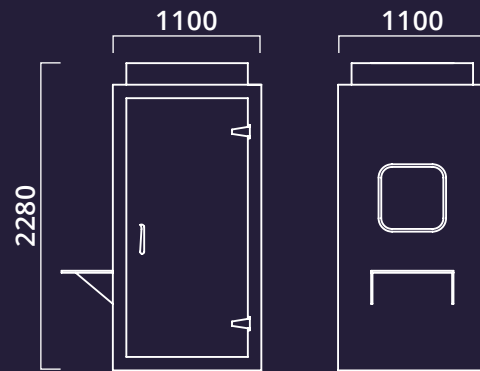

Superior noise reduction


- Extra thick 100mm panels for high noise reduction
- Double-glazed window, quality acoustic seals, double-wall steel construction
- Wall cavities packed with high density, sound-absorbent acoustic padding
- Designed to meet Australian Standard AS1269.4 specification
- Factory tested for internal noise levels prior to dispatch
- Perforated steel interior reduces reverberation (echo)

Quality construction

- Solid construction using quality Australian steel
- High-quality components, fittings, seals and paint
- Manufactured to suit a range of locations (mine sites, medical centres, research labs, etc.)
- Custom-designed, heavy duty door hinges
- Designed and built by Australian acoustic fabricators with over 30 years' experience

Spacious and durable

- Largest booth in the single person range
- Extra leg and elbow room to accommodate people of all shapes and sizes
- Easy entry/exit: low step, wide door opening, high ceiling
- Very stable, floor-mounted structure
- Add window in door for extra light and visibility

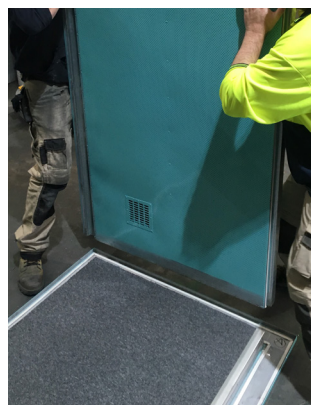

Tongue and groove assembly



Mains powered light and fan assembly

Design features

- Flat-packed tongue and groove design, for straightforward onsite assembly
- Can be disassembled, relocated and reassembled
- Standard 6.35mm sockets suit most audiometers
- External shelf (add or remove as required)
- Shelf slot for easy connection and less cable clutter
- Choose left or right door opening to suit room layout
- Wide selection of colours in hard-wearing finish
- Wired to Australian Standards specification with circuit breaker and safety switch

Which booth is right for you?

	Ezi Booth (most popular)	Flexi Booth	Maxi Booth
External dimensions	Width: 1050mm Depth: 810mm Height: 2000mm	Width: 1050mm Depth: 980mm Height: 2100mm	Width: 1100mm Depth: 1100mm Height: 2280mm
Internal dimensions	Width: 890mm Depth: 640mm Height: 1670mm	Width: 950mm Depth: 840mm Height: 1910mm	Width: 900mm Depth: 900mm Height: 1940mm
Wall thickness	50mm	50mm	100mm
Weight	350kg	250kg	500kg
Assembly	Delivered fully assembled	3 hours, 2-3 people	4 hours, 3-4 people
Portable	On castors - can be wheeled between rooms or buildings	Disassembly required for relocation	Disassembly required for relocation

Power	240 volt 10 amp (mains power), circuit breaker, safety switch
Door opening	Choose left or right opening
Inclusions	5m extension cable, external shelf, internal fan, low-heat LED light, carpet, equipment hooks
Options	Extended warranty, maintenance kit, USB socket, additional window in door

For sales and service in your area:

Glenn Johnson

☎ 1300 118 222

✉ contact@hearingcompany.com.au